POINT-OF-CARE

osom® COVID-19/Flu A&B Combo Test

OSOM® COVID-19/Flu A&B Combo Test

ONE SAMPLE AND TEST TO SIMULTANEOUSLY DIFFERENTIATE BETWEEN COVID-19, FLU A, AND FLU B

Features and Benefits

- **Reliable:** Proven clinical performance to help you reduce the risk of misdiagnosis & inappropriate treatment.
- **Rapid:** Differentiate between three viral infections (COVID-19, Flu A, & Flu B) in only 15 minutes, so you can facilitate immediate treatment protocols at the first patient visit.
- **Cost Efficient:** Unique QC Inside® feature which provides you with two additional tests per kit for external QC testing.
- **Simple Workflow:** An easy CLIA-waived procedure that can be performed in only a few steps.
- **Superior Supply Chain:** Manufactured in the USA so you can avoid delays from international product shipments.

Performance Characteristics

	Sensitivity*	Specificity*	Sample Type	Limit of Detection**	Storage Temperature	Time to Result
SARS-CoV-2	87.6%	99.8%	Anterior nasal swab	1.27×10^4	15-30°C/ 59-86°F	15 minutes
Flu A	86.9%	99.6%		1.4×10^5		
Flu B	86.8%	99.7%		2.25×10^5		

*Combined prospective clinical RT-PCR comparator study.

**LOD for SARS-CoV-2 BA.5 Omicron Variant. Refer to IFU for additional LODs.




QC inside

	FLU A	FLU B	SARS-CoV-2
SENSITIVITY	86.9%	86.8%	87.6%**
SPECIFICITY	99.6%	99.7%	99.8%


SEKISUI DIAGNOSTICS

Because every result matters™

OSOM® COVID-19/Flu A&B Combo Test

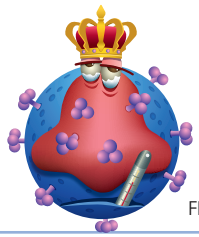


Three results in only 15 minutes

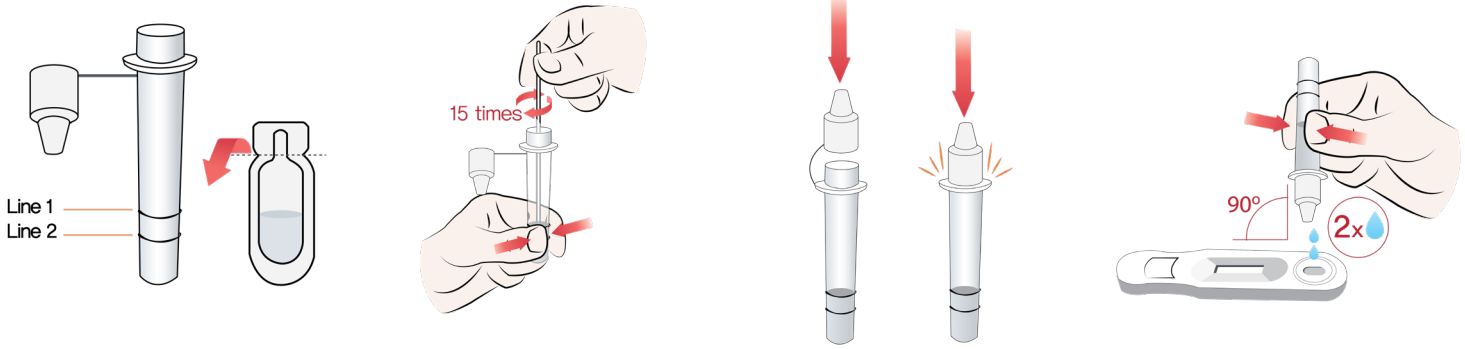

1. **FILL** the sample tube with reagent solution above line 1.
2. **INSERT** the swab into the sample tube and swirl at least 15 times.
3. **CLOSE** the affixed dropper tip of the sample tube firmly.
4. **ADD** 2 drops of the sample to the sample well of the test device and wait 15 minutes for results.

Kit includes one Flu A/Flu B/SARS-CoV-2 positive Control Swab and one Flu A/Flu B/SARS-CoV-2 negative Control Swab

INTERPRETATION OF RESULTS

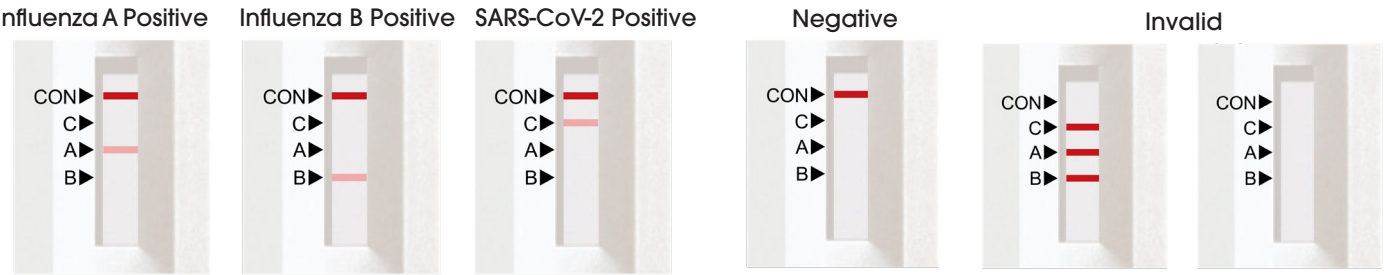

Ordering Information

	Configuration	Catalog Number	CPT Code
OSOM COVID-19/FLU A&B COMBO TEST	25 tests per kit	1099	87812QW

To learn more about all of our OSOM Tests scan the QR code or go to osomtests.com

SEKISUI DIAGNOSTICS
Because every result matters™
sekisuidiagnostics.com

MADE IN THE USA

© 2026 SEKISUI Diagnostics, LLC. All rights reserved. OSOM® is a registered trademark of SEKISUI Diagnostics, LLC.
Because every result matters™ is a trademark of SEKISUI Diagnostics, LLC.

80-8977-00-03 09/26